



Options and Accessories Visuals

Base

Flat Free Tires - Block and V-Profile



V-Profile for drive and front wheel for indoor use



Block-Profile for drive and front wheel for outdoor use

V-Track Headrest (incl. Track extension)



V-Trak Headrest winged
 V-Trak Headrest simple

V-Track Thoracic support & pad



V-Trak Thoracic support & pad & cover large right & left
 V-Trak Thoracic support & pad & cover medium right & left

V-Track Backrest/Standard & Axxis



Offers more back and neck support

Control Options

R-Net Joystick Module



R-Net includes connector points for accessory "Buddy Button"

R-Net Omni+ Module



Compatible with many input devices
 Remote control
 4/5 switch control
 3 Axis switch or proportional head control
 Sip & Puff
 Single switch scanner
 (For use with Flexible Arms or Joint Arm)

Easy Touch



Alternative keypress control with dual function positioning

Swing away joystick holder



Allows for easier transfer
 Provides opportunity to maneuver chair in close proximity of tables

Joystick protection guard



Protects hands & Joystick

LCEV Stand-up mechanism



Standard on standing operation (armrest round with control switch)

Control Options

LCEV Armrest



Armrest 6 cm/2.4" wide with built-in control switch
 Armrest 6 cm/2.4" wide with hand control
 Armrest round with hand control
 Armrest 5 cm/2" longer

Swing away chest support



Rotatable and swing away to the rear. Increased feeling of safety for the user in standing position. Distance between backrest tube and the beginning radius

LCEV Dual switch for table control



Stand-up operation with push buttons for table. Table tray not included

Buddy Button



Green: ON/OFF
 Yellow: MODE
 Red: Emergency stop

LCEV Control switch extension



Vertical for round armrest
 Vertical for wide armrest
 Horizontal

Joystick-top



Customized control grips

T-bar Mushroom
 Stick Chin
 Softball Ball

Seating System

Powered reclining back/leg rest 90-180°



Powers to complete prone/recline

Powered recline and stand



Powers from lying position directly to standing

Tilt in Space



Tilts 0-35 degrees

Leg Rest Options

Knee wedge

Rectangular



Knee supports are only delivered extra strong. Stabilization is removable

Conical



This knee wedges give more support, especially with spasticity

Leg Rest Options

Knee support - PRO




The new Knee support "PRO" option provides enhanced stability/rigidity and allows for separate positioning of each knee/leg

Knee support - Integral



Individually adjustable

One-piece footrest fix, no flip-up feature

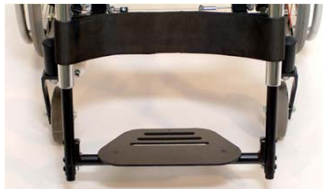

Option- Easier to handle more stability

Ground level footrest



Especially for users who use the wheelchair as a stand-up help

Two-piece footrest



Standard- Easier to move in and out of chair

Footplate protection frame



Protection frame only. Protects toes from collision injuries

Footrest height adapter



Especially for users with short lower leg length. The footrest becomes 1.2" narrower than the seat width

Seat depth / Footrest quick adjustment



This feature makes adjustment on footplate height quick and easy

Suitable for demo chairs

Accessories

Curb Climber



Front castor fork with double wheel. For curb climbing up to 3.2" / 8.1 cm

Foam seat cushion with cover



Seat cushion option

Backrest Extension - 6" (15CM)



Only available for standard Tension Adjust Sling Back backrest

Back rest extension with push handle and upholstery



Length:
5.9"
9.8"
13.8"

Accessories

Stroller handle/Push Handle Bar



Can be used with standard & V-Track back. Helps for manual pushing of the chair

Table tray



Fixed tube in armrest

Swing away table tray



Support installed at the armrest. Table is still fastened on the chair, even when not in use
Only for Combi arm rest

Swing away table tray center cut joystick holder



Depth adjustable. Table is still fastened on the chair, even if not in use. Joystick can be hinged away, if not used

Skirt guard



Easy to remove for lateral transfer. (Pair)

Thigh support



Supports the thigh and stabilizes the leg axle in seating and standing position

Chest role



Individually adjustable in height. Offers additional support but does not replace the chest strap

Foot guide



Can be adjusted individually. Secures and aligns feet

Footrest wedge



Offers foot alignment for cases of spasticity and special need

Chin control with shoulder holder



Shoulder- and chest holder with upholstery. Height & width are individually adjustable using "ON/OFF" & "MODE" button

Chin control swing away SPZ



Electrically operated. Height, width & angle adjustable

Attendant module



Attendant controls
Easy operation transition

Arms to adapt the R-Net Omni+ Module



Flexible arm including adapter
Joint arm including adapter stronger and more adjustable
Can be bent in every direction

Monitor



Provides an objective outcome measurement for users and therapists
Delivered with PC-Software.

Accessories

Rear wheels with drum breaks



Allows assistant to stop the wheelchair, safe and comfortably. Cannot be combined with the option of Spinergy, Camber or One arm drive

Rear wheels "Spinergy"



Operation without additional energy. Black, blue, red, yellow & white

Steel hand ring



Increased grip, less temperature-sensitive and less dirty hands

Maxgrepp (Hand ring with rubber covered)



Original

Better hold for persons with limited strength or grip



Ergo Para CP

One arm drive



For stroke patients or people with unilateral limited force or grip
Only possible with standard wheel

Push handle



Standard push handles

Push handle extension adjustable



Helpful if the wheelchair is pushed by different persons

Traditional Anti-tip rolls



Anti-tip rolls to keep chair from tipping backwards. Right & Left sold separately

Anti-tip aid / tip aid swing away "HC"



Anti-tip aid to keep chair from tipping backwards which includes easy to use swing away function. Right & Left sold separately

Brake lever extension (6 cm/2.4")



Simpler accessibility and less strength necessary

Rear wheels with hub locking brakes



Light weight, improved. Cannot be combined with Hub locking brake, Drum brakes or One arm drive

Spoke guard



Spoke guard-Transparent

Spoke guard-Transparent w/LEVO-Logo

Vehicle Mount & Driving Accessories

Vehicle tie-down kit



Secure tie-down for vehicle transport

EZ-Lock brackets



Levo chairs adjust to EZ-Lock

Lights and indicators (traffic approved)



Turn signals, brake lights, & lamps for on road use

Rear view mirror



Mounted on Joystick side or mounted opposite side

Mudguard for drive wheel



Mud Guards/Sold in pairs. Silver; Black & Rubber

Frame reinforcement for e-fix



Wheelchairs with E-Motion need frame reinforcement. You can only have this combined with anti-tip aid "HC" and not with anti-tip rolls or tip aid