

Etac Clean Mobile Shower Commode

 etac®

Clean is a prize winning best seller because of its many smart functions.



Clean – the Smart Shower Chair

Improved Access

The seat has a recess at the front and is fully open at the rear. This improves the accessibility for intimate hygiene.

The foot support is always available but never in the way.

The foot support has a unique design. The support can slide in underneath the seat when not needed. Thereby it is not lost and always at hand when requested. This is also handy for transfers and the armrests which can be taken off or swung away will be of additional help.



Clean fits both floor-mounted and wall-mounted toilets.

Adjustable between 47.5–60 cm (18¾"–23½") 2.5 cm (1") increments.



Easy Maneuvering

Clean is very easy to maneuver because it's dainty and the shape of the handles which makes the chair easy to turn. There is also a strong advantage in the fact that Clean fits both floor and wall mounted toilets perfectly.

Durable Construction

Clean has neither welding seams nor holes which reduces the risk of corrosion.

Clean – Height Adjustable

The height can easily be adjusted without any tools to fit perfectly over the toilet, or to simplify transfers.

Clean 24" is easy to Maneuver


Clean 24" is a well-balanced shower chair that allows the user independence. The seat is fully open at the rear and facilitates for the user to manage intimate hygiene unassisted.

Side transfers can be performed safely thanks to sufficient seat space in front of the wheels. Furthermore the brake handles are high mounted without obstructing transfer.


With a rear wheel adapter, the regular Clean with a 55 cm (21¾") seat height can be changed into a Clean 24". There are many accessories to choose from also for Clean 24".




Clean

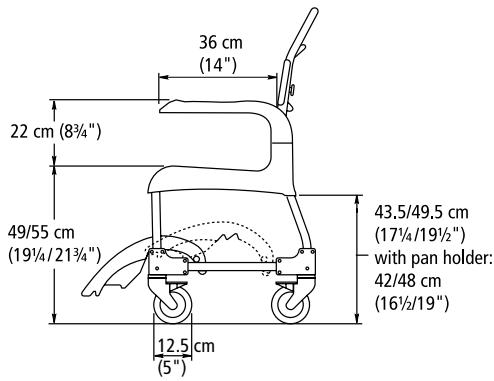

Clean height adjustable

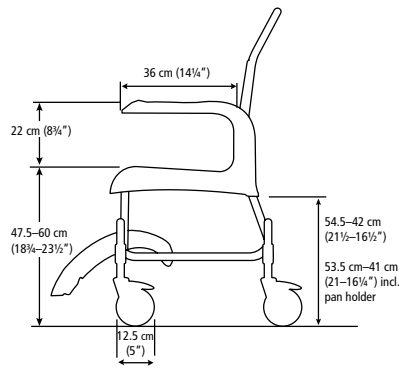


Clean 24" rear wheels

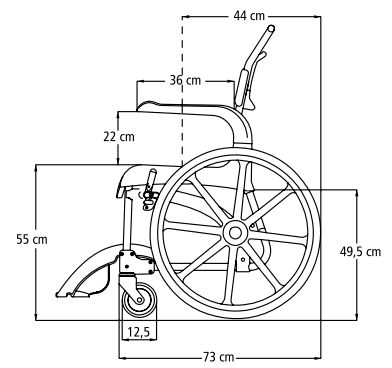

Other dimensions

	Clean	Clean height adjustable	Clean 24" rear wheels
Total width	52 cm (20½")	55 cm (21¾")	69 cm (27½")
Width between armrests	43.5 cm (17¼")	43.5 cm (17¼")	43,5 cm (17¼")
Seat width	48 cm (19")	48 cm (19")	48 cm (19")

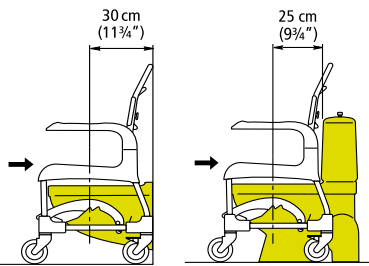

Clean is suitable for both floor and wall mounted toilets.

Description

	Item no. Green	Grey
Clean 44 cm (17¼") 4 lockable castors	8020 9209	8020 9269
Clean 49 cm (19¼") 4 lockable castors	8020 9210	8020 9270
2 lockable castors	8020 9213	8020 9273
4 lockable castors, double box	8020 9220	---
2 lockable castors, double box	8020 9223	---
Clean 55 cm (21¾") 4 lockable castors	8020 9211	8020 9271
2 lockable castors	8020 9214	8020 9274
4 lockable castors, double box	8020 9221	---
2 lockable castors, double box	8020 9224	---
Clean 60 cm (23½") 2 lockable castors incl heightening spacers	8020 9215	9020 9275
Clean height adjustable	---	8020 9208
Clean 55 cm (21¾") with 24" rear wheels	8020 9216	8020 9276

Materials

Chassis: Powder-painted steel tubing
 Seat/backrest: Polypropylene
 Armrests: Polypropylene
 Footplate: Polypropylene
 Castor brackets: Polyamide
 Castors, lockable, ball-bearing: Polyamide/TPE/Rust free tubing
 Rear wheels: Polyamide, mounted in ball bearings

Weight

Clean: 14 kg (30.8 lbs)
 Clean height adjustable: 14,2 kg (31.3 lbs)
 Clean with 24" rear wheels: 20 kg (44.1 lbs)

Care

Regular cleaning plus inspection of castors

Design

A&E Design



Clean the product with a solvent-free cleaning agent (pH level 6-9), or with a 70% disinfectant solution. Can be decontaminated at max. 85°C / 185°F if so required.



Accessories

Pan

Green 8020 9230
Grey 8020 9254



Lid with handle

Green 8020 9232
Grey 8020 9256

Pan incl. lid with handle

Green 8020 9231
Grey 8020 9255

Material: Polyamide

Pan holder

Green 8020 9233
Grey 8020 9257
Fits both Clean pans and regular hospital pans with a thin rim.



Bucket holder

Green 8020 9234
Grey 8020 9258
Fits bucket with a thick rim.

Material: Polypropene

Attachment for low mounting of pan/bucket

Green 8020 9235
Grey 8020 9259
Fitted to the pan/bucket holder



Material: Polypropene

Soft pad

Green 8020 9236
Grey 8020 9260
Heat-reflecting soft seat with good anti-slip properties even when wet.
Recess width: 17 cm (6¾")



Material: Cellular polyethylene

Cushion, Protect 8020 9245

Air-inflated RoHo cushion with excellent pressure-relieving properties. The air pressure is adjustable.
Recess width: 21 cm (8¼")
Material: Neoprene



Comfort pad 8020 9248

A comfortable, hygienic seat pad which makes the hole a little smaller. Easy to keep clean.
Recess width: 18 cm (7")



Material: Polyurethane

Seat cover 8020 9247

When a complete seat is required. Easy to keep clean.



Material: Polyurethane



Soft back pad

Green 8020 9238
Grey 8020 9261

Heat-reflecting soft back pad which provides extra comfort for the user.
Material: Cellular polyethylene



Soap dish

Green 8170 1140
Grey 8020 9266

Attached to the side of the chair. Keeps the soap within convenient reach.
Built-in holder for the shower handle.



Cross bar 8020 9246

Attached to the armrests.
Material: Aluminium, cellular polyethylene, polypropylene



Safety belt 8020 9034

Attached to the rear of the backrest.
Material: Polyester



Heightening spacers 8020 9243

Raise the seat by five centimetres. Do not fit Clean with 24" rear wheels or Clean height adjustable.
Material: Aluminium

80 mm lockable castors, 2 pcs.

Green 8020 9240
Grey 8020 9263



80 mm non lockable castors, 2 pcs.

Green 8020 9241
Grey 8020 9264

Reduce seat height by 5 cm (2"). Mounted on Clean 44 cm as standard. Do not fit Clean with 24" rear wheels or Clean height adjustable.

Material: polyamide/TPE



Rear wheel adapter 8020 9242

Only fits Clean with a seat height of 55 cm (21¾").
Material: Aluminium, polyamide



Angle adjustable leg rest

Item no. (right) 8020 9252
Item no. (left) 8020 9253

Attachments included. Only fits Clean height adjustable.
Material leg support: ABS plastic



Folding footrest, pair

Item no. 8020 9250

Removable, adjustable height. Attachments included. Only fits Clean height adjustable.
Material footplate: Polypropylene